# Propaq® Encore Monitor

---

---

## Update: Directions for Use/Reference Guide

English	Update to Directions for Use/Reference Guide: Software Version 2.5X
Deutsch	Aktualisierungen/Ergänzungen zur Bedienungsanleitung: Softwareversion 2.5X
Español	Actualización de la Guía de referencia: Versión de software 2.5X
Français	Mise à jour du Guide de référence : Logiciel Version 2.5X
Italiano	Aggiornamento delle Guida di riferimento: versione 2.5X
Nederlands	Update voor de Gebruiksaanwijzing: softwareversie 2.5X
Português	Actualização das Instruções de Utilização : Versão de Software 2.5X
Svenska	Uppdatering av bruksanvisning: Programversion 2.5X

**WelchAllyn®**

Advancing Frontline Care™

Copyright 2005, Welch Allyn. All rights are reserved. No one is permitted to reproduce or duplicate, in any form, this manual or any part thereof without permission from Welch Allyn.

Welch Allyn assumes no responsibility for injury or for any illegal or improper use of the product that may result from failure to use this product in accordance with the instructions, cautions, warnings, or indications for use published in this manual.

Welch Allyn®, Propaq® Encore, Smartcuf®, Acuity®, FlexNet® and Flexible Monitoring® are registered trademarks of Welch Allyn, Inc.

Nellcor® is a registered trademark of Nellcor Puritan Bennett. Masimo® and SET® are registered trademarks of Masimo Corporation.

Software in this product is Copyright 2005, Welch Allyn or its vendors. All rights are reserved. The software is protected by United States of America copyright laws and international treaty provisions applicable worldwide. Under such laws, the licensee is entitled to use the copy of the software incorporated with this instrument as intended in the operation of the product in which it is embedded. The software may not be copied, decompiled, reverse-engineered, disassembled or otherwise reduced to human-perceivable form. This is not a sale of the software or any copy of the software; all rights, title and ownership of the software remains with Welch Allyn or its vendors.

For information about any Welch Allyn product, call Welch Allyn Technical Support:

USA	800 535 6663 315 685 4560	Australia	61 2 9638 3000 800 074 793
Canada	800 561 8797	China	86 21 6327 9631
European Call Center	353 46 9067790	France	33 1 55 69 58 49
Germany	49 7477 9271 86	Japan	81 3 3219 0071
Latin America	305 669 9003	Netherlands	31 157 505 000
Singapore	65 6419 8100	South Africa	27 11 777 7555
United Kingdom	44 1332 363812		

### Reorder Number 810-1719-XX

Manual Part Number 810-2107-00 Rev. A 6/2005

For a hard-copy version of the *English Propaq Encore Reference Guide* and the English section of this Update, order P/N 810-2106-XX.

Welch Allyn  
8500 Creekside Place  
Beaverton, Oregon 97008-7107  
USA  
[www.monitoring.welchallyn.com](http://www.monitoring.welchallyn.com)

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# English Update: Propaq Encore Monitor

---

This document describes changes in operation for Propaq Encore monitors (software version 2.5X). This document replaces information in the following manuals:

- *Propaq Encore Reference Guide* (software version 2.4X)
- *Propaq Encore Directions for Use* (software version 2.4X)

For monitor operating information not contained in this document, please refer to the preceding manuals.

This document also describes basic information about how Propaq Encore monitors operate when they are connected to an Acuity Central Monitoring System. For additional information about the Acuity System, refer to the Directions for Use for your Acuity System.

## Intended Use

*The following statement is an addition to the Intended Use section in the Propaq Encore Directions for Use and the Propaq Encore Reference Guide.*

Propaq Encore monitors that do not include CO<sub>2</sub> or printer options are able to withstand light rain exposure over short periods of time (uniform distribution of approximately 1 mm of water/ minute for 10 minutes or less).

# Symbols

## Symbols in This Update

**WARNING** Indicates conditions or practices that could result in illness, injury or death

**Caution** In this manual, indicates conditions or practices that could damage the equipment or other property.

**Note** Provides additional information.

## Symbols on the Propaq Encore Monitor, Accessories and Packaging

	Non-ionizing electromagnetic radiation		<b>Caution:</b> On the product, means "Consult accompanying documentation."
	Separate lead acid battery from other disposables for recycling		Fuse
Pb			
	This device has been tested and certified by the Canadian Standards Association International to comply with applicable U.S. and Canadian medical safety standards.		Single use only
	Temperature limits		Stacking limit (by number)
	Altitude limit		Humidity limit
	Keep away from rain		Fragile

# Warnings

**WARNING** When taking NIBP measurements, periodically observe the patient's limb to make sure that the circulation is not impaired for a prolonged period of time. Also make sure the blood pressure cuff is properly placed according to the *Propaq Encore Directions for Use* or *Propaq Encore Reference Guide*. Be especially careful when using the short-term automatic mode (TURBOCUFF). Prolonged impairment of circulation or improper cuff placement can cause contusions.

**WARNING** The range of values measured by the monitoring parameters is provided in the Specifications section of the *Propaq Encore Directions for Use* or *Propaq Encore Reference Guide*. Operation of the monitor outside the range of specified values is not recommended and may cause inaccurate results.

**WARNING** Electronic equipment that emits very strong electromagnetic or radio frequency signals can cause electrical interference with monitor operation, including causing the monitor to turn off power. Avoid operating this monitor near such equipment. For guidance about electromagnetic emissions and the recommended separation distance between the monitor and such equipment, refer to the specifications section of this manual.

**WARNING** For best product performance and measurement accuracy, use only accessories supplied by Welch Allyn or recommended in the Welch Allyn *Products and Accessories* booklet (P/N 810-0409-XX). Use accessories according to your facility's standards and the manufacturer's recommendations. Always refer to the manufacturer's Directions for Use. Use of accessories not recommended by Welch Allyn may result in increased electromagnetic emissions or decreased electromagnetic immunity of the equipment.

# Patient Alarms and Equipment Alerts

## Description of Alarm and Alert Tone Patterns

### A New Tone for Patient Alarms

- Patient alarm tone pattern: one second on, two seconds off

### Other Alarm and Alert Tones

- Apnea alarm tone pattern: one second on, one second off (fastest)
- Equipment alert tone pattern: one second on, four seconds off (slowest)

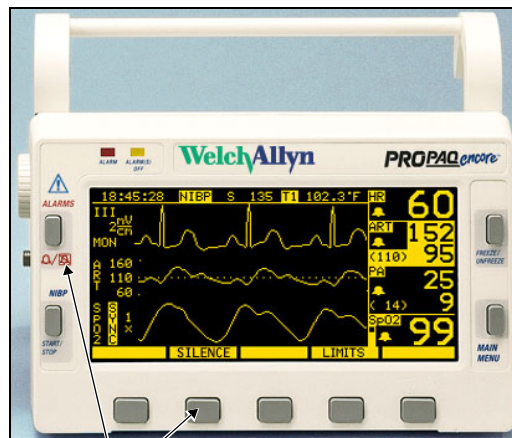
**Note** Propaq monitors connected to an Acuity Central Monitoring System can sound distinct tone patterns for certain alarms, such as life-threatening Arrhythmia and ST alarms. For further information regarding Acuity System alarm tones, please consult the *Acuity System Directions for Use*.

## Silence an Active Patient Alarm or Equipment Alert Tone for 90 Seconds

You can silence the tone of a patient alarm or equipment alert for a period of 90 seconds.

### Silence an Alarm or Alert Tone for 90 Seconds

Press either key to silence the alarm or alert tone for 90 seconds.

1. Check the patient and provide appropriate care.
2. Press or the **SILENCE** key.
3. After caring for the patient, make sure the appropriate alarm limits are set.

### Re-enable an Alarm or Alert Tone Before the 90-Second Silence Period Has Elapsed

1. Press .

### Indications During a 90-Second Silence Period

At the Propaq Encore Monitor	At the Acuity Central Station
During the silence period, the usual <i>visual</i> alarm and alert indications are displayed.	
If a new alarm or alert occurs, the new alarm or alert tone interrupts the silence period.	Only life-threatening arrhythmia and apnea alarm tones interrupt the suspend period.
To re-enable an alarm or alert tone before the 90-seconds has elapsed, press  .	To resume an alarm or alert tone before the 90-seconds has elapsed, press <b>RESUME</b> .
If an equipment alert occurs during the silence period, you can acknowledge (dismiss) <i>all</i> indications of the alert by pressing any key, as usual.	
If the original alarm or alert condition still exists after 90 seconds, the alarm or alert tone sounds again.	

## Inhibit Alarm and Alert Tones for Four Minutes: 4 SUSPND

As you perform patient care, there might be occasions when you want to suspend potential or current alarm and alert tones for a period of four minutes.

### Initiate a Four Minute Suspend Period

1. Press **MAIN MENU, SETUP, ALARMS, 4 SUSPND**.

### Resume Alarm and Alert Tone Capability Before the Four Minute Period Has Elapsed

1. Press /  .

### Indications During a Four-Minute Suspend Period

At the Propaq Encore Monitor	At the Acuity Central Station
If an alarm or alert occurs during the suspend period, the usual <i>visual</i> alarm and alert indications are displayed.	
<i>Propaq monitors connected to Acuity Central Station:</i> Apnea or life-threatening arrhythmia alarm tones interrupt the suspend period.	Apnea or life-threatening arrhythmia alarm tones interrupt the suspend period.
<i>Standalone Propaq monitors:</i> Apnea alarms <i>do not</i> interrupt the suspend period, and arrhythmias are not detected.	
To resume alarm and alert tone capability before the 4 minutes has elapsed, press  /  .	To resume alarm and alert tone capability before the 4 minutes has elapsed, press <b>RESUME</b> .
If an equipment alert occurs during the suspend period, you can acknowledge (dismiss) <i>all</i> indications of the alert by pressing any key, as usual.	

## Inhibit Alarm and Alert Tones Indefinitely: ALL ALARMS

**Note** This feature requires access to the Service Menu. The Service Menu features are not intended for use during ordinary, routine operation.

You can suspend potential or current alarm and alert tones for an indefinite period of time. The tones are disabled until someone re-enables them.

### Indefinitely Suspend All Alarm and Alert Tones

1. Press **MAIN MENU, SETUP, MORE, MORE, SERVICE, YES, MORE, MORE, SETTINGS, ALL ALARM, YES.**

**WARNING** Whenever audible alarm tones are disabled, make sure the patient is closely observed.

### Resume All Alarm and Alert Tone Capability

1. Press / .

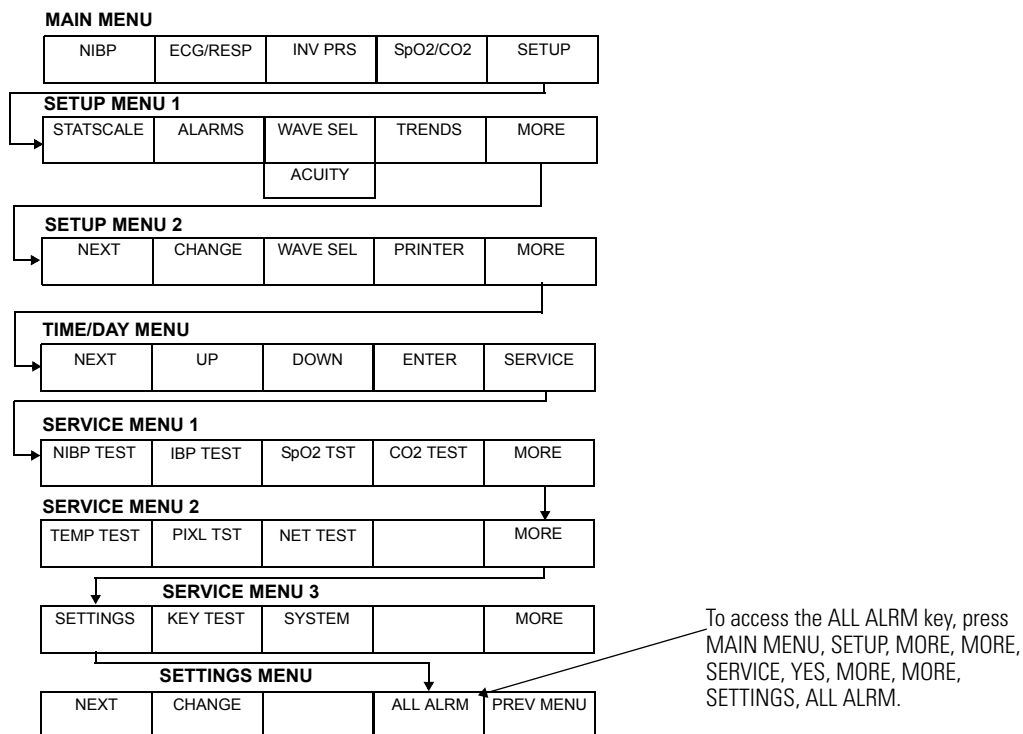
### Indications During an All Alarms Suspended Period

At the Propaq Encore Monitor	At the Acuity Central Station
If an alarm or alert occurs during the suspend period, the usual <i>visual</i> alarm and alert indications are displayed.	
<i>Propaq monitors connected to Acuity Central Station:</i> Apnea or life-threatening arrhythmia alarm tones interrupt the suspend period.	Apnea or life-threatening arrhythmia alarm tones interrupt the suspend period.
<i>Standalone Propaq monitors:</i> Apnea alarms <i>do not</i> interrupt the suspend period, and arrhythmias are not detected.	
To resume alarm and alert tone capability, press  /  .	To resume alarm and alert tone capability, press <b>RESUME</b> .
If an equipment alert occurs during the suspend period, you can acknowledge (dismiss) <i>all</i> indications of the alert by pressing any key, as usual.	

## New Location for the ALL ALRM Key

The ALL ALRM key has been relocated to the Settings Menu under the Service Menu.

## Respond to Equipment Alerts

**Note** If multiple equipment alerts occur simultaneously, the Propaq monitor displays a MULTIPLE EQUIPMENT ALERT message. In certain situations, such as when a channel is disabled, the source of the second alert is not displayed in the new message. Look for blank numeric and waveforms areas on the screen to identify the source(s) of multiple alerts.

## Summary of Alarm and Alert Keys and Acuity Central Station Messages

The following table summarizes Silence and Suspend behaviors and Acuity Central Station messages for different combinations of current and previous Propaq monitor and Acuity software versions.

**Note** When your facility uses both older and newer Propaq monitors and Acuity software versions, be aware that messages vary at Acuity Central Station. Some Acuity Central Station messages do not identify specific alarm/alert tone silence or suspension states.

However, visual alarm indications still continue as usual at the Acuity Central Station and at the Propaq monitors.

Propaq Monitor Version	Propaq Monitor Key	Result of Propaq Monitor Key Press	At Acuity Central Station, Result of Propaq Monitor Key Press	
			Acuity 6.30 and prior	Acuity 6.31 and later
Propaq Encore 2.5X	SILENCE or 	Silences active alarm and alert tones for 90 seconds at the Propaq monitor and Acuity Central Station.  All new alarms and alerts interrupt the silence.	Alarms Suspended Message  Only life-threatening arrhythmia and apnea alarms interrupt the silence.	Alarms Suspended Message  Only life-threatening arrhythmia and apnea alarms interrupt the silence.
	4 SUSPND	Inhibits alarm and alert tones for four minutes at the Propaq monitor and Acuity Central Station. <sup>a</sup>	Alarms Suspended Message	Nurse Suspend 4 min Message
	ALL ALRM	Indefinitely inhibits alarm and alert tones at the Propaq monitor and Acuity Central Station, until someone resumes tone capability. <sup>a</sup>	Alarms Suspended Message	Nurse Suspend Always Message
Propaq Encore 2.4X and prior	SUSPEND or 	Suspends (for 90 seconds) the capability for patient alarm and equipment alert tones at the Propaq monitor and Acuity Central Station. <sup>a</sup>	Alarms Suspended Message	Alarms Suspended Message
	ALL ALRM	Turns off all patient vital sign alarm limits at the Propaq monitor and Acuity Central Station.  No visual or auditory alarm indications occur until someone resumes the limits. <sup>a</sup>	Some Alarms Off, Set Alarm Limits Message	Some Alarms Off, Set Alarm Limits Message

a. If connected to the Acuity Central Station, life-threatening arrhythmias and apnea alarms interrupt the suspension.

## Confirm and Learn Alarm Behavior in the Inservice Mode

To confirm that a Propaq monitor is properly generating patient alarms and to learn the alarm behavior of the Propaq Encore monitor, perform the following steps.

1. Disconnect all patient cables from the Propaq monitor.
2. To temporarily remove customized alarm settings, press the **MAIN MENU** key, then **SETUP, MORE, CHANGE, SETUP, USE NOW, YES**.
3. To put the Propaq monitor in (non-alarming) Inservice Mode 1, press **MAIN MENU, SETUP, WAVE SEL, INSERV**.
4. To set the Propaq monitor in (alarming) Inservice Mode 2, press **MAIN MENU, SETUP, WAVE SEL, INSERV**.

Within 3-5 seconds, the monitor will alarm because the “patient’s” vital signs fall outside of the alarm limit ranges.

5. Confirm the following:

The monitor is sounding an alarm  
A vital sign numeric is flashing  
The small red light in the top right corner of the monitor is flashing  
The SILENCE and LIMITS keys are showing on the screen

6. To silence the alarm tone for 90 seconds, press /  or **SILENCE**.

Visual alarm indications remain, and the alarm tone resets after the 90-second silence period.

7. To bring the Propaq monitor back into the (non-alarming) Inservice Mode 1, press **INSERV**.

The “patient’s” vital sign readings will return to acceptable levels within the alarm limit range.

8. Turn off the Propaq monitor.

Any previously set custom patient mode settings are restored when you turn it on again.

Under these conditions, if the monitor fails to generate visual or auditory alarm indications, carefully repeat the above steps. If the monitor is still unresponsive, remove it from circulation and take it to your facility’s biomedical service department for evaluation.

## Avoid Electrostatic Discharge

When humidity in the working environment decreases, the human body and other insulators can become charged with static electricity due to friction.

To prevent unwanted electrostatic discharge (ESD), follow these standard guidelines:

- Maintain the recommended humidity of 40% to 60% in the work environment.
- Dissipate electrostatic charge before performing routine operator maintenance.

## Specifications

### Monitor (Physical) Specifications

Characteristic	Specification
Recovery time following defibrillator discharge	Less than or equal to 10 seconds

### Monitor (Environmental) Specifications

Characteristic	Specification
Degree of protection against ingress for monitors <i>without</i> CO <sub>2</sub> or printer options	IPX1 rating, drip proof per EN60529: 1991

## Electrosurgery

Propaq Encore monitors are suitable for use in the presence of electrosurgery.

## EMC Compliance

Special precautions concerning electromagnetic compatibility (EMC) must be taken for all medical electrical equipment. The Propaq Encore complies with IEC EN 60601-1-2:2001.

- All medical electrical equipment must be installed and put into service in accordance with the EMC information provided in this document and the *Propaq Encore Directions for Use* or *Propaq Encore Reference Guide*.
- Portable and mobile RF communications equipment can affect the behavior of medical electrical equipment.

Propaq Encore monitors comply with all applicable and required standards for electromagnetic interference.

- It does not normally affect nearby equipment and devices.
- It is not normally affected by nearby equipment and devices.
- It is safe to operate the monitor in the presence of high-frequency surgical equipment.
- However, it is good practice to avoid using the monitor in extremely close proximity to other equipment.

---

### Propaq Encore Monitor: Guidance and Manufacturer's Declaration—Electromagnetic Emissions

---

The Propaq Encore monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the monitor should assure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic Environment—Guidance
RF emissions CISPR 11	Group 1	The monitor uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The monitor is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

---

**Propaq Encore Monitor: Guidance and Manufacturer's Declaration—Electromagnetic Immunity**

The Propaq Encore monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the monitor should assure that it is used in such an environment.

<b>Immunity Test</b>	<b>IEC 60601 Test Level</b>	<b>Compliance Level</b>	<b>Electromagnetic Environment—Guidance</b>
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations on power-supply input lines IEC 61000-4-11	<5% $U_t$ (>95% dip in $U_t$ ) for 0.5 cycle 40% $U_t$ (60% dip in $U_t$ ) for 5 cycles 70% $U_t$ (30% dip in $U_t$ ) for 25 cycles <5% $U_t$ (>95% dip in $U_t$ ) for 5 sec	<5% $U_t$ (>95% dip in $U_t$ ) for 0.5 cycle 40% $U_t$ (60% dip in $U_t$ ) for 5 cycles 70% $U_t$ (30% dip in $U_t$ ) for 25 cycles <5% $U_t$ (>95% dip in $U_t$ ) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the monitor requires continued operation during power mains interruption, it is recommended that the monitor be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note	$U_t$ is the AC mains voltage prior to application of the test level.		

**Propaq Encore Monitor: Guidance and Manufacturer's Declaration—Electromagnetic Immunity**

The Propaq Encore monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the monitor should assure that it is used in such an environment.

<b>Immunity Test</b>	<b>IEC 60601 Test Level</b>	<b>Compliance Level</b>	<b>Electromagnetic Environment—Guidance</b>
			Portable and mobile RF communications equipment should be used no closer to any part of the monitor, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
			<b>Recommended Separation Distance</b>
Conducted RF IEC 61000-4-6	3 V <sub>rms</sub> 150 kHz to 80 MHz	3 V <sub>rms</sub>	$d = 1.2 \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	$d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz <p>where <b>P</b> is the maximum output power rating of the transmitter in watts according to the transmitter manufacturer and <b>d</b> is the recommended separation distance in meters.  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey<sup>a</sup>, should be less than the compliance level in each frequency range<sup>b</sup>.  Interference might occur in the vicinity of equipment marked with the following symbol:</p> 
Note 1	At 80 MHz and 800 MHz, the higher frequency range applies.		
Note 2	These guidelines might not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.		
a	Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the monitor is used exceeds the applicable RF compliance level above, the monitor should be observed to verify normal operation. If abnormal performance is observed, additional measures might be necessary, such as reorienting or relocating the monitor.		
b	Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.		

---

**Recommended Separation Distances Between Portable and Mobile RF Communications Equipment and the Propaq Encore Monitor**

---

The Propaq Encore monitor is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the monitor can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the monitor as recommended below, according to the maximum output power of the communications equipment.

Rated Maximum Output Power of Transmitter W	Separation Distance According to Frequency of Transmitter m		
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance  $d$  in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where  $P$  is the maximum output power rating of the transmitter in watts (w) according to the transmitter manufacturer.

---

Note 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

---

Note 2 These guidelines might not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

---

## Patient Alarm and Equipment Alert Specifications

Characteristic	Specification
System alarm delay time	0.5 seconds, typical
Indicators	Flashing red ALARM light: patient alarm(s) Continuously on red ALARM light: patient alarm silenced or suspended  Flashing yellow ALARMS OFF light: equipment alert Continuously on yellow ALARMS OFF light: one or more patient alarm limit(s) off
Tone patterns	Apnea: 1 second on, 1 second off Patient alarm: 1 second on, 2 seconds off Equipment alert: 1 second on, 4 seconds off
SILENCE and  /  keys	Silences patient alarm or equipment alert at the Propaq monitor and Acuity Central Station for 90 seconds or until  /  is pressed again.  New patient alarms or equipment alerts interrupt the silence at the Propaq monitor, but only life-threatening arrhythmia and apnea alarm tones interrupt the silence at Acuity Central Station.
4 SUSPND key	At the Propaq monitor and Acuity Central Station, inhibits all patient alarm and equipment alert tone capability for 4 minutes or until  /  or RESUME is pressed.  Apnea and life-threatening arrhythmia alarm tones do interrupt the silence for Propaq monitors connected to the Acuity Central Station.
ALL ALRM key	At the Propaq monitor and Acuity Central Station, inhibits all patient alarm and equipment alert tone capability indefinitely, until  /  or RESUME is pressed.  Apnea and life-threatening arrhythmia alarm tones do interrupt the silence for Propaq monitors connected to the Acuity Central Station.
RESUME and  /  keys	At the Propaq monitor and Acuity Central Station, re-enables capability for alarm and alert tones and re-enables alarm and alert tones that are still active.
Acknowledge equipment alert	Tone and visual message are dismissed.
HR/PR alarm limits range	25 to 250 beats per minute
Apnea delay alarm limit range	Neonate: 6 to 20 seconds Adult/Pediatric: 6 to 30 seconds
IBP alarm limit ranges	All patient modes Systolic, Diastolic, Mean -30 to 300 mmHg
NIBP alarm limit ranges	Neonate Systolic 25 to 120 mmHg Diastolic 10 to 105 mmHg Mean 10 to 110 mmHg Pediatric Systolic 30 to 160 mmHg Diastolic 15 to 130 mmHg Mean 15 to 140 mmHg Adult Systolic 30 to 260 mmHg Diastolic 20 to 235 mmHg Mean 20 to 255 mmHg
Temp alarm limit range (T1, T2)	32.0 to 122.0 °F 0.0 to 50.0 °C
Δ Temp alarm limit range	0.0 to 90.0 °F 0.0 to 50.0 °C

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